



## Using your Bluetooth Mobile Phone with the Accent 700 Quick Reference Guide

Although there are many phones to choose from these days, the phone you use depends on whether you would like to just make and receive voice calls or also be able to send and receive text messages. If purchasing a new phone, we recommend you go into the phone retailer and test to ensure that the phone is compatible. We have found **Sony Ericsson** phones to be most compatible.

### PAIRING YOUR MOBILE PHONE

To pair your phone with your device you need to make sure your **Mobile Phone** has **Bluetooth** switched **ON** and that the phone is **Discoverable**. To connect the phone to the device:

1. Go to **Toolbox**.
2. Select **Maintenance Menu**.
3. Select **Hardware Diagnostics**.
4. Select **Phone Configuration**.
5. Select **Pair To Phone**.
6. A yellow box will appear on the device advising "**Please Wait**", and then a few seconds later a message will appear on the phone asking to "**Add to devices**", select **Yes**.
7. Enter the pass code "**0000**" when prompted.
8. The Bluetooth Status Box will display "**Phone Connected & Text Messaging Supported.**"



Once **paired**, the **PHONE CONNECTION** icon appears in the **Status Display Area**. The icon goes **green** and displays the **battery charge** and **signal strength** of the **Mobile Phone**.



9. Select **OK** twice to exit the **Phone Configuration Menu**.

### PHONE ACTIVITY

For quick access to your phone functions, a **Phone Activity** has been created for you and is present in the **2<sup>nd</sup> Row of your Activity Row**.



### MAKING A TELEPHONE CALL FROM THE DIAL PAD


The **Phone Dial Pad** allows you to: **Dial calls, Answer calls, use your Core vocabulary to converse, Send and Receive Text Messages**.



1. Select **Phone Dial Pad**. Type in the number in. \*\* When dialing a number you may or may not need to add "**0044** or **44**" in front of the number, for e.g. **0044** or **4479810662517**. This is solely to do with how networks register the number via Bluetooth.\*\*
2. Select the **Send** button
3. **GO TO CORE**, to enable you to converse.
4. You may also go to the **PHONE PAGE**, where you will find some conversational phrases, or to **SPELL** if you wish to use **SPELLING**.
5. To end the call, go back to the **Phone Dial Pad** and select **End**.

**Note:** Use comma's when you dial phone numbers that contain extensions. For example: **3302621984,123**. Use multiple commas for multiple extensions.

### ANSWERING A CALL FROM THE DIAL PAD

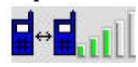1. Upon receipt of a call, the **CALL STATUS** icon appears  in the **Status Display Area** of your device.
2. Select the **Send** button on the **Phone Dial Pad** to accept the call.
3. Select the **Go to Core** key and use your **Core Vocabulary** to converse.
4. To end the call, go to the **Phone Dial Pad**, and select **End**.


\*\* The **Phone** activity contains a **Phone Send/End** key which allows you to answer and hang-up without having to use the **Dial Pad**.



### USING THE INTERNAL MICROPHONE

There is an **Internal Microphone** built in the device that you may use to speak on the phone if you are able to do so. The option is located within the **Phone Configuration Menu**.

1. Go to the **Phone Dial Pad**, and select **Phone Configuration**.
2. Turn the **Microphone** option to **ON**.
3. Adjust the **Mic Level** option to make the volume louder or softer. \*\*If the volume is too loud, you will hear feedback noise.\*\*
4. You may go to **CORE** to converse whether the microphone is on or not.

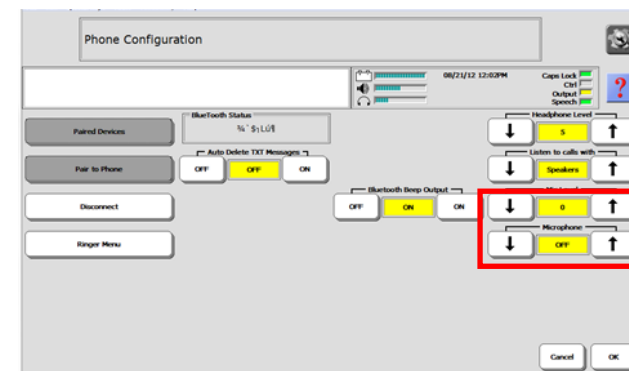

## SENDING TEXT MESSAGES

1. Select **Send Text Message** from the **Phone Dial Pad**. Or **PHONE SEND TEXT MESSAGE** from within the **Phone Activity**.
2. Enter the telephone number you wish to send the message to. **\*\*Remember Do not use Spaces or Hyphens.** You may need to add **"0044 or 44"** in front of the number.\*\*
3. Select **OK**.
4. You will be told to **"Enter Text Message,"** type your text message and select **OK**.
5. The device will notify you with the following, **"Sending Message"** and then **"Message Sent Successfully."** Select **OK**.

## SENDING MESSAGES TO EXISTING CONTACTS

1. Select the **Phone Dial Pad** from within your **Activity Row**.
2. Select **Contacts** and choose the Contact you wish to send a **Text Message** to.
3. You will be taken back to the **Phone Dial Pad**. Select **Phone Send Text Messages**.
4. You will be told to **"Enter Text Message,"** type your text message and select **OK**.
5. The device will notify you with the following, **"Sending Message"** and then **"Message Sent Successfully."** Select **OK**.

## READING TEXT MESSAGES

Upon receipt of a Text Message, a **white envelope icon** appears in the **Status Display Area**.

1. Select **Phone Text Messages** from within the **Phone Activity**.
2. The device will display a list of text messages from your phone. Select the most recent.
3. Read the message, and then select the desired action at the bottom of the screen: **Reply, Delete, Speak Message, Add Contact, Copy Message to Clipboard** or **Cancel**

## ADDING CONTACT NAMES & TELEPHONE NUMBERS TO THE PHONE DIAL PAD

1. In the **Phone Dial Pad**, select **Add Contact**.
2. Type the phone number you wish to add. **Do not use Spaces or Hyphens.** \*\* You may need to add **"0044 or 44"** in front of the number, for e.g. **004479810662517**.

3. Select **OK**.
4. Type the **Name** of the **Contact**. Select **OK**.
5. This takes you back to the **Phone Dial Pad**.
6. Select **Contacts**, to view the added **Contact**.


## ADDING AN INCOMING NUMBER


When you receive a call, the caller's number will be displayed in the **Text Display Area**. If you haven't already, you can add this number to your **Contacts**.

1. Go to the **Phone Dial Pad**.
2. Select the **Add Contact** button. A keyboard with the **Number** displayed will appear.
3. Select **OK**.
4. Enter the **Name** of your **Contact**.
5. Select **OK**. You will be taken back to the **Phone Dial Pad**.
6. Select **Contacts**, to view the added **Contact**.

## SAVING CONTACTS TO PHONE DIRECTORY PAGE

Once you have your **Contacts** in the **Phone Dial Pad**, you may wish to add these to your **Phone Directory Page** - accessed via the **Phone Activity**, this will allow you to have your very own personalized **Phone book**.

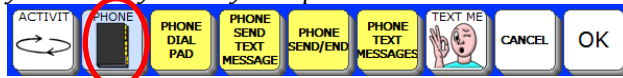

1. Select **PHONE** from within the **Phone Activity**.
2. Press the tools key  then select the **SETUP KEY** from the on-screen user keys
3. Select one of the **Name** keys. You may also select a **Blank** key if you wish.
4. Select **Spell Message or Define Key Function**
5. Press **Clear Display**.
6. Select **Insert Tool** and scroll through by pressing **More Items** to select **Phone Dial Number**.
7. Your **Phone Dial Pad Contacts** will be displayed. Select the **Contact** you wish to add.
8. Select **OK**.
9. Select **Change Icon** – you may wish to use an imported image of the **Contact**.
10. Change the **Label** of the icon to your **Contact's Name**.
11. Repeat steps 2 to 10 to add more contacts.

## SAVING YOUR CONTACTS LIST

**Important Note:** To save your **CONTACTS**, you **MUST** perform an **Entire Device Contents** save and not a **One User Area** save – for saving data from the **Accent 1200** see the **Accent 1200's general Quick Reference Guide**. Alternatively, you could save your contacts separately.

1. Insert your **USB** stick into the device.
2. Go to **Toolbox**.
3. Select **Transfer Memory Menu**.
4. Select **Save Contacts**. The **"Save"** window will appear on the screen.
5. Type a name for your save. The file type, **.mti** will be displayed. **Do not delete the file Type**.
6. Select **OK**.
7. A yellow box tells you when the transfer is complete. Select **OK**.
8. Select **OK** to exit.

## LOADING YOUR CONTACTS

1. Insert your **USB** stick into the device.
2. Go to **Toolbox**.
3. Select **Transfer Memory Menu**.
4. Select **Load Contacts**. \*\*A yellow box will appear advising that the current contacts will be erased, if you have added contacts since your last save, you may not want to load an earlier list. In that case select **NO** in the yellow box. If you want to continue, select **YES**. \*\*
5. Select the **.mti** file you wish to load from your **USB** stick.
6. Select **OK** twice to exit.

## Liberator contacts

### Regional Consultants

Mark Street: 07747016660

Eleri Davies: 07827324286

Jacqui Malthouse 07825512669

### Office

Technical Support: 01733 370470

(Option 2)